


















JOVIVE WALNUT CREEK - WOUND CARE DRESSING GUIDE








Quick reference for wound care products available in clinic

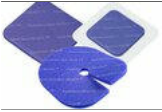


PRIMARY DRESSING OINTMENTS AND GELS			
Image	Name & Type	Description & Instructions	Frequency
	Medihoney Gel <i>Medical-grade honey gel</i>	Promotes debridement and reduction of bioburden. Available on unit carts.	Change: Daily
	Hydrogel <i>Dry wounds</i>	Adds moisture to dry wounds. Available on unit carts.	Change: Daily
	Iodosorb <i>Cadexomer iodine gel</i>	Absorbs drainage and promotes debridement. Comes from pharmacy.	Change: 2-3 times weekly
	Betadine <i>Short-term infection control</i>	Broad-spectrum antiseptic. Promotes eschar formation, prevents MASD. Apply sparingly. Can delay healing if used excessively.	Change: As needed
	Mupirocin <i>Localized bacterial skin infections</i>	Topical antibiotic effective against gram-positive bacteria.	Change: 2-3 times daily
	Gentamicin <i>Infected wounds with gram-negative organisms</i>	Topical antibiotic effective against gram-negative bacteria.	Change: 2-3 times daily
	Zinc Oxide <i>Peri wound protection, moisture barrier</i>	Protective barrier cream for moisture-prone areas.	Change: As needed
	Nystatin Topical Powder <i>Fungal infections, moisture areas</i>	Antifungal powder for treating candida infections in skin folds.	Change: 2-3 times daily










	Lidocaine 5% Ointment <i>Topical anesthesia for painful wounds</i>	Local anesthetic ointment for pain relief during wound care.	Change: As needed
Each of these can be spread over a piece of oil emulsion gauze first, then applied to the wound, to ensure coverage of the wound bed.			
PERI-WOUND TREATMENTS			
Image	Name & Type	Description & Instructions	Frequency
	Triad Cream <i>Incontinence-associated dermatitis</i>	Protects periwound and intact skin from moisture damage. Ideal for sacral breakdown. Reapply after each soiling.	Change: After each soiling
	Calmoseptine <i>Minor skin irritation</i>	Small tube format of moisture barrier and antiseptic ointment.	Change: As needed
	Antifungal Cream <i>Fungal skin infections</i>	Topical antifungal for treating superficial fungal infections.	Change: 2-3 times daily
	Hydrocortisone Cream 2.5% <i>Inflammatory skin conditions</i>	Mid-potency topical steroid for inflammatory dermatoses.	Change: 1-2 times daily
	Triamcinolone Acetonide 0.1% Cream <i>Inflammatory dermatoses</i>	Mid-potency topical steroid for eczema and inflammatory skin conditions.	Change: 1-2 times daily
	Aquaphor <i>Dry skin, friction areas, minor wounds</i>	Occlusive healing ointment for dry and damaged skin. Criticaid replacement.	Change: As needed
	Gold Bond <i>Dry skin, irritation protection</i>	Medicated powder and cream for skin protection and moisture. Sween replacement.	Change: As needed
	Amlactin <i>Dry, rough skin</i>	Lactic acid moisturizing lotion for dry and rough skin. Urea replacement.	Change: Daily









PRIMARY DRESSINGS










Image	Name & Type	Description & Instructions	Frequency
	Hydrofera Blue <i>Methylene Blue/Gentian Violet Foam</i>	Blue foam, absorbs drainage. Plain blue side touches wound. Obtain extra from SWAT.	Change: 1-3 days
	Adaptic/Telfa Oil Emulsions <i>Fragile skin, graft sites, skin tears</i>	Nonadherent contact layer. Use when minimizing trauma during dressing changes is important. Apply over ointments if needed.	Change: Daily
	Fibracol Collagen Dressing <i>Chronic wounds, granulation stimulation</i>	Collagen matrix dressing that supports wound healing and tissue regeneration.	Change: Weekly
	Xeroform <i>Petroleum gauze with Bismuth</i>	Non-adherent, antimicrobial sheets. Useful for skin tears, burns, painful wounds. On unit carts. Replacement: oil emulsion.	Change: Daily
	Calcium Alginate (Plain) <i>Moderate to heavy exudate</i>	Non-silver alginate dressing for absorbing wound drainage.	Change: 1-3 days
	Calcium Alginate with Silver <i>Infected wounds with heavy exudate</i>	Silver-impregnated calcium alginate for antimicrobial wound management.	Change: 1-3 days
	Gelling Fiber with Silver (Rope) <i>Deep tunneling wounds with infection risk</i>	High-absorbent silver gelling fiber in rope form for tunnels and undermining.	Change: 1-3 days










SECONDARY DRESSINGS




Image	Name & Type	Description & Instructions	Frequency
	Non-adherent <i>Moist, superficial wounds</i>	Nonstick surface for moist, superficial wounds. Often paired with secondary gauze wrap.	Change: Daily

	ABD Pads <i>Heavily draining wounds</i>	Layer of absorption. Not a primary dressing. Used over packing or as pressure padding.	Change: As needed
	Kerlix Wrap <i>Securing dressings to limbs/head</i>	Stretch gauze; allows breathable coverage. Can be used with saline as wet-to-dry dressing base.	Change: As needed
	Foam/Silicone Bordered <i>Moderate drainage wounds</i>	Cushions & absorbs moderate drainage, water resistant outside.	Change: 2-3 days
	Super Absorbent Dressing <i>Heavily draining wounds</i>	Absorbs excess amounts of drainage/blood.	Change: 1-2 days
	Spandagrip Compression Size E <i>Light compression support</i>	Graduated compression bandage for support and edema control.	Change: Daily
	Spandagrip Compression Size F <i>Moderate compression support</i>	Medium compression bandage for moderate edema and support.	Change: Daily
	Spandagrip Compression Size G <i>Firm compression support</i>	Higher compression bandage for significant edema management.	Change: Daily
	Hydrocellular Non-Adherent Foam <i>Moderate to heavy exudate</i>	Non-adherent foam dressing that requires secondary dressing.	Change: 1-3 days
	Thin Hydrocolloid Dressing <i>Shallow wounds, minimal exudate</i>	Thin hydrocolloid for superficial wounds with light drainage.	Change: 3-7 days
	Retention Tape 2 inch <i>Securing large dressings</i>	Non-woven retention tape for securing bulky or large dressings.	Change: As needed

	Cohesive Bandage (Tan) <i>Self-adhering wrap</i>	Self-adhering elastic wrap that sticks to itself, not skin.	Change: As needed
	Conform Stretch Bandage <i>Securing dressings</i>	Lightweight stretch bandage for securing dressings without constriction.	Change: As needed
	Island Border Gauze Dressing <i>General wound coverage</i>	Sterile bordered gauze dressing with adhesive perimeter.	Change: Daily
	Tubular Elastic Net Size 6 <i>Dressing retention on limbs</i>	Tubular elastic net for holding dressings in place on arms/legs.	Change: As needed
	Tubular Elastic Net Size 9 <i>Dressing retention on torso</i>	Larger tubular elastic net for torso and larger body areas.	Change: As needed
MULTI-LAYER COMPRESSION WRAPS			
Image	Name & Type	Description & Instructions	Frequency
	Urgo K2 Dual Compression System <i>Ankle compression 30-40mmHg</i>	Two-layer compression system for venous insufficiency and edema management.	Change: Weekly
	Urgo K2 Lite Dual Compression <i>Light compression therapy</i>	Lighter compression system for sensitive patients or mild edema.	Change: Weekly
	Unna Boot with Calamine <i>Venous leg ulcers, edema control</i>	Zinc oxide-impregnated bandage with calamine for venous insufficiency.	Change: Weekly
<i>Used for venous insufficiency and edema control. Ensure arterial circulation is adequate before application.</i>			
ADDITIONAL SUPPLIES AND TOOLS			
Additional Wound Care Supplies			

	Saline Irrigation Solution <i>Wound irrigation and cleaning</i>	Sterile 0.9% sodium chloride for wound irrigation and cleaning.	Change: As needed
	Dermal Wound Cleanser Spray <i>Wound irrigation and cleaning</i>	No-rinse wound cleanser spray for gentle wound bed preparation.	Change: As needed
	Vashe Wound Wash <i>Wound cleaning and irrigation</i>	Hypochlorous acid wound cleanser that's gentle yet effective.	Change: As needed
Instruments & Tools			
	Silver Nitrate Applicators <i>Cauterization, granulation control</i>	Caustic applicators for cauterizing small bleeding points and hypergranulation.	Change: As needed
	Disposable Biopsy Punch <i>Tissue sampling</i>	Sterile disposable punches for skin biopsy procedures.	Change: Single use
	Dermal Curette 2mm <i>Tissue debridement</i>	Disposable curette for gentle tissue debridement and sampling.	Change: Single use
	Dermal Curette 4mm <i>Tissue debridement</i>	Disposable curette for tissue debridement of larger areas.	Change: Single use
	Surgical Blade #15 Carbon Steel <i>Precise tissue cutting</i>	Sterile disposable surgical blade for debridement and minor procedures.	Change: Single use
	Skin Staple Remover Tray <i>Staple removal</i>	Sterile staple remover with prep materials for safe staple removal.	Change: Single use

	Suture Removal Tray <i>Suture removal</i>	Sterile tray with instruments for safe suture removal procedures.	Change: Single use
	Paper Wound Measuring Ruler <i>Wound measurement documentation</i>	Disposable paper rulers for accurate wound measurement and documentation.	Change: Single use
Prep & Cleaning Supplies			
	Adhesive Remover Wipes <i>Safe adhesive removal</i>	Gentle removal of adhesive residue from skin without trauma.	Change: As needed
	Alcohol Prep Pads <i>Skin antisepsis</i>	70% isopropyl alcohol for skin preparation and antisepsis.	Change: As needed
	Povidone Iodine Prep Pads <i>Antiseptic skin preparation</i>	10% povidone iodine prep pads for surgical site preparation.	Change: As needed
	Benzoin Tincture Compound <i>Adhesive enhancement</i>	Enhances tape/dressing adhesion and protects skin from moisture.	Change: As needed
	No-Sting Skin Prep Wipes <i>Skin protection before adhesive application</i>	Protective barrier wipes that enhance adhesion without stinging.	Change: As needed
	Sanicloth Super Germicidal Wipes <i>Surface disinfection</i>	Broad-spectrum disinfectant wipes for equipment and surface cleaning.	Change: As needed
Basic Supplies			
	Gauze Sponge 2x2 (8-ply) <i>Wound cleaning, packing</i>	Small sterile gauze sponges for wound care and cleaning.	Change: Single use

	Gauze Sponge 4x4 (12-ply) <i>Wound cleaning, absorption</i>	Standard sterile gauze sponges for wound care and cleaning procedures.	Change: Single use
	Cotton-Tip Applicators 6 inch <i>Precise application, wound cleaning</i>	Sterile cotton-tipped applicators for precise medication application.	Change: Single use
	Tongue Depressor Blades <i>Oral examination, mixing</i>	Sterile wooden tongue depressors for examinations and mixing compounds.	Change: Single use
	Reinforced Skin Closure Strips <i>Wound edge approximation</i>	Sterile adhesive strips for closing minor lacerations and incisions.	Change: Until healing
	Paper Adhesive Tape 2 inch <i>Gentle dressing securement</i>	Hypoallergenic paper tape for sensitive skin applications.	Change: As needed
Specialty Items			
	Iodoform Nasal Packing Strip <i>Nasal hemorrhage control</i>	Antiseptic-impregnated packing strips for nasal bleeding control.	Change: Single use
	Plain Nasal Packing Strip <i>Nasal support, non-antiseptic packing</i>	Non-medicated packing strips for nasal support and hemostasis.	Change: Single use
	Brava Moldable Ostomy Ring <i>Ostomy seal enhancement</i>	Moldable ring for improving ostomy appliance seal and preventing leakage.	Change: With appliance changes
	Saline Inhalation Solution <i>Respiratory therapy</i>	0.9% saline for nebulization and respiratory treatments.	Change: As prescribed
General Supplies			

Important: Always assess wound bed, surrounding skin, and patient factors before selecting dressing. For detailed protocols visit: <https://github.com/bobburgers7/wound-care-dressings>

Change Log

Version 0.3.7 - 2025-09-08

- Added Xeroform back to PRIMARY DRESSINGS section
- Moved from removed items to active product list
- Positioned between Fibracol Collagen and Calcium Alginate Plain

Version 0.3.6 - 2025-09-04

- Reorganized wound care product sections:
 - Renamed "COMMONLY ORDERED OINTMENTS AND GELS" to "PRIMARY DRESSING OINTMENTS AND GELS"
 - Renamed "COMMONLY ORDERED WOUND FILLERS AND MEDICATED DRESSINGS" to "PRIMARY DRESSINGS"
 - Renamed "COMMONLY ORDERED COVER DRESSINGS" to "SECONDARY DRESSINGS"
 - Created new "PERI-WOUND TREATMENTS" section
 - Moved Triad, Calmoseptine, antifungal cream, hydrocortisone cream, and triamcinolone to peri-wound section
- Removed generic "Calmoseptine" entry and renamed "Calmoseptine Ointment 1/8oz" to "Calmoseptine"
- Removed "collagen" from primary dressing ointments section
- Added three new peri-wound treatment products:
 - Aquaphor (Criticaid replacement)
 - Gold Bond (Sween replacement)
 - Amlactin (Urea replacement)
- Created optimized images:
 - gold-bond.jpg (59% size reduction: 162KB → 66KB)
 - amlactin.jpg (68% size reduction: 74KB → 24KB, converted from .jpeg)
- Fixed page break issue in product grid: section headers now stay with their products when printing

Version 0.3.5 - 2025-08-29

- Fixed font warnings by removing unsupported "sans-serif" fallback from font declarations
- Improved changelog spacing in PDF output for better readability
- Optimized changelog formatting in typ files so bullet list entries weren't so far apart

Version 0.3.4 - 2025-08-29

- Updated PDF filenames to include ISO date format (YYYY-MM-DD) at the beginning
- Updated QR code and reference links to point to GitHub repository: <https://github.com/bobburgers7/wound-care-dressings>
- Added date consistency requirements to CLAUDE.md instructions

Version 0.3.3 - 2025-08-29

- Renamed product names for clarity:
 - "Iodasorb / Betadine" → "Betadine"
 - "Telfa Pads" → "Non-adherent"
- Created optimized zinc-oxide.jpg image (72% size reduction)

Version 0.3.2 - 2025-08-29

- Added changelog display to both main.typ and product-grid.typ documents
- Updated product-grid.typ to use configuration-based date and version (matching main.typ)
- Fixed compression image references to use correct .jpg extensions
- Created optimized compression-k2.jpg and compression-k2-lite.jpg images
- Added Pacific Time Zone notation to CLAUDE.md instructions
- Updated document version synchronization between config and changelog

Version 0.3.1 - 2025-08-28

- Removed multiple wound care products from active guide and moved to removed-items file:
 - Ointments/Gels: Santyl, Silvadene, Vaseline, Aquaphor, Wound Hydrogel (Amorphous)
 - Filler/Medicated: Medihoney Alginate, Gelling Fiber, Xeroform (noted oil emulsion as replacement), Prisma Promogran, Silvasorb, Aquacel Ag, Alginates/Silver, Hydrofera Blue ReadyTransfer, Xtrasorb Super Absorbent
 - Cover Dressings: Allevyn, Tegaderm Absorbent, DuoDerm, ABD & Rolled Gauze, Coban/ACE, Netting/Stockinette, Tegaderm Film, Bordered Gauze/Island Dressing, Hypafix Retention Tape, Allevyn Life Lite 4x4, Allevyn Life Lite Large, Silicone Foam with Border 6x6, Silicone Foam with Border 3x3
- Updated configuration files to reflect product removals
- Maintained all removed items in separate JSON file for reference

Version 0.3.0 - 2025-08-23

- Major refactoring: unified single JSON data source, implemented smart image optimization, added multi-layer compression wrap category, reorganized additional supplies with subcategories, improved image alignment in grid layout, removed unnecessary catalog items, created automated version management system

Version 0.2.0 - 2025-08-23

- Added multi-layer compression wrap category, unified single JSON data source, optimized images for print PDFs, added organized additional supplies section